

endolite
get busy living

aqualimb

Unique fast assembly
system for shower limbs

The water resistant Aqualimb prosthesis is designed
for safe use in showers and wet rooms.



www.endolite.co.uk



AqualimbTT:



Pre-formed cosmetic shin for ease of finishing



Max. user weight:	80kg system 100kg system
Activity level:	1-2
Size range:	22cm-27cm
Component weight:	1.39kg* 100kg
Build height:	280 to 415mm 100kg 100 to 415mm 80kg
Heel height:	5mm
Fitting instruction:	80kg - 938123 100kg - 938255
Warranty:	24 months


*Component weight shown is for a size 24cm

100kg Part number:

Size	Left	Right
22cm	629705L	629706L
23cm	629707L	629708L
24cm	629709L	629710L
25cm	629711L	629712L
26cm	629713L	629714L
27cm	629715L	629716L

80kg Part number:

Size	Left	Right
22cm	629405LF	629406LF
23cm	629407LF	629408LF
24cm	629409LF	629410LF
25cm	629411LF	629412LF
26cm	629413LF	629414LF
27cm	629415LF	629416LF

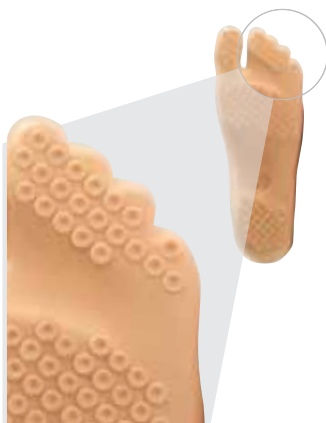

The Aqua Foot anti-slip tread pattern is designed for safe barefoot walking on wet surfaces

The Aqualimb is designed for use by all amputees for showering and bathing purposes, it has an Activity Level 2 rating only.



Example

AQTF 25L
Size/Side

Foot example: aqualimbTF, size 25 left

Max. user weight:	100kg
Activity level:	2
Size range:	22cm-27cm
Component weight:	2.0kg
Build height:	400-585mm
Heel height:	5mm
Instructions for use:	938321
User guide:	938326
Warranty:	24 months

*Component weight shown is for a size 24cm

AqualimbTF:



The Aqua Knee is a moulded polycentric linkage design for improved stability. It also has a manual lock option for increased security when needed.

The Fast assembly package has a cut and slide connection for length adjustment and domed Endolite style alignment for gait optimisation.



endolite
get busy living

